

Antigen-Down HRP Conjugate Stabilizer, 5X


Preserves the activity of HRP conjugates and reduces background signal.

Antigen-Down HRP Conjugate Stabilizer, 5X is supplied as a 5X concentrate. When diluted to its 1X working concentration, it can be used to create concentrated stock conjugates for long-term storage, reconstitute lyophilized HRP conjugates, and dilute IgG-peroxidase conjugates to their useful working titer in ELISAs and other immunology-based techniques.


At 1X, Antigen-Down HRP Conjugate Stabilizer preserves the functional integrity of the horseradish peroxidase enzyme and immunoglobulin components of the HRP-IgG conjugate complex. It utilizes a hydrolyzed, non-mammalian protein matrix instead of the conventional BSA protein stabilizers used in other commercial stabilizing products. Incorporation of a heterogeneous, smaller molecular weight protein matrix may allow for creation of more intimate associations between stabilizer proteins and the performance-critical regions of the IgG-HRP conjugate, resulting in vast improvements in stabilizing performance. In addition, the physical and chemical characteristics of the protein matrix create an unfavorable environment for non-specific conjugate binding interactions with the ELISA plate surfaces. When assessing low-positive, humoral antibody responses in an antigen-down ELISA format, such as to a food allergen or early signs of a viral pathogen outbreak in a swine herd, minimization of non-specific background signal is of critical importance and can easily be the difference between success or failure in the assay development process.

By maintaining conjugate component activity and preventing contamination issues through the inclusion of an antimicrobial agent, Antigen-Down HRP Conjugate Stabilizer extends the functional utility of stored and reconstituted IgG-peroxidase conjugates. Additional components may be added to address unique ELISA sample situations or for use in other immunological techniques.

Stabilizes HRP Conjugates

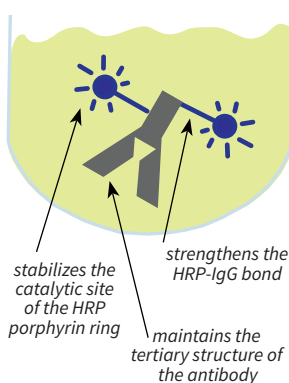

Minimizes Background in Antigen-Down ELISAs

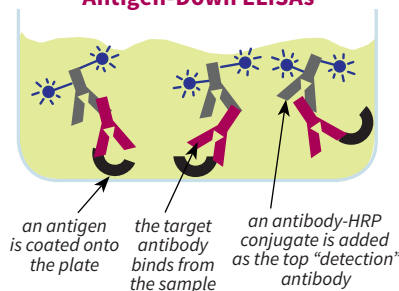

ANTIGEN-DOWN HRP CONJUGATE STABILIZER, 5X

Size	Catalog#
25 mL	LS-M30-25
100 mL	LS-M30-100
500 mL	LS-M30-500
1 L	LS-M30-1

INSTRUCTIONS:

1. Protect the stabilizer from light.
2. Gently mix the stabilizer; avoid bubbles.
3. Prepare 1X working strength stabilizer solution by diluting 1:5 in diH₂O. Add 1 part 5X stabilizer to 4 parts diH₂O.
4. Mix for 15 minutes.
5. Prepare the HRP conjugate solution at the predetermined dilution factor or concentration.
6. Mix for 15 minutes.
7. Pipette 50-300 μ L of the diluted HRP conjugate to each well of the ELISA at the appropriate step in the assay.
8. Store remaining conjugate solution at 2-8°C. Protect from direct exposure to light.

For more ELISA protocols and information, please visit www.LSBio.com.

SPECIFICATIONS:

- Clear to light yellow liquid
- 5X concentrate
- pH 7.2-7.6

STORAGE:

- 24 months at 2-8°C
- 1 week at room temperature

SAFETY & USAGE:

- Warning! May cause an allergic skin reaction. Harmful to aquatic life with long lasting effects.
- SDS available upon request.
- Product intended for research use or for further manufacturing into *in vitro* diagnostics reagents only.
- Not intended for use in human or therapeutics purposes.



2401 4th Avenue, Suite 900, Seattle, WA 98121

USA TOLL-FREE 1-866-819-4732 LOCAL 206-374-1102

www.LSBio.com